

SereneLife

Questions? Issues?
We are here to help!
Phone: (1) 718-535-1800
Email: support@pyleusa.com

SereneLife



PSLMSGBMB22 Thermal Therapy Massage Kit

Innovative Waterless Heating Technology
with Bamboo Sticks, and Heating Bag





User Manual


PLEASE MAKE SURE THAT YOU HAVE READ THE MANUAL CAREFULLY BEFORE USING IT.

THERMAL THERAPY MASSAGE KIT

Thermal therapy is the use of heat in therapy, such as for pain relief and health. It can take the form of a hot stone, hot bamboo sticks, and others. It can be beneficial to those with arthritis and stiff muscles and injuries to the deep tissue of the skin. Hugeworth born out of the desire to make massage therapy mobile. Our innovative electric massage warmer set comes with a light-weight energy efficient stylish storage/travel bag case that doubles as stone heater oven that quickly heats up the rocks, and it is the helpful tool to heat the massage stones or bamboo sticks. As we know the application of warmed Basalt stones/bamboo sticks combined with manual massage to provide the client with a deeply relaxing Therapeutic experience. This treatment which also benefits the therapist not only financially but also physically.

PRODUCT FEATURE


-  **Innovative waterless heating technology**
With waterless heating technology there is no need for traditional stone heating equipment. It can drop the electricity consumption, much "greener" and more environmentally friendly.
-  **Digital temperature controller**
Our digital temperature controller grants therapist the ability to choose the desired stone temperature.
-  **Portable massage kit**
Our portable massage kit offers mess-free stone warming along with the ability to easily transport for traveling therapist as well as your massage practice or home use.
-  **Multifunctionality**
Combined heating with bag, it is more than a bag, not only can self-contained but also can heat.

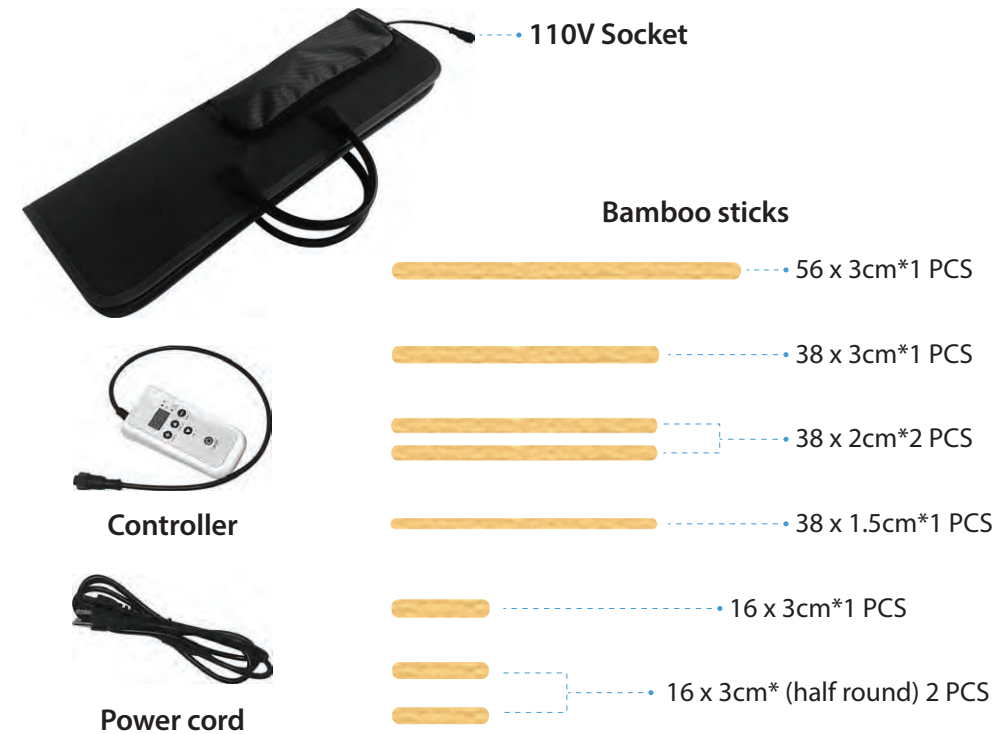

HOW TO USE


STEP 1

Place Bamboo sticks inside the portable massage heating warmer. Make sure all are placed flat and not stacked on top of one another.

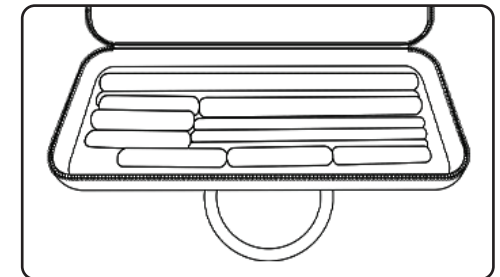


STEP 2

Connect the temperature controller to the Thermal Therapy Massage Kit.

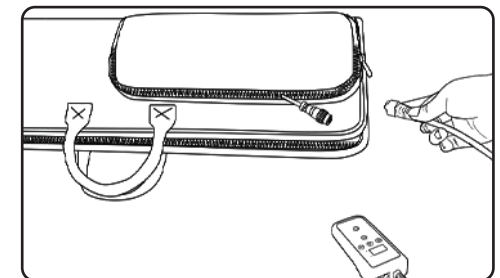


STEP 3

Connect the power cord to the temperature controller and then plug power cord into a power outlet.

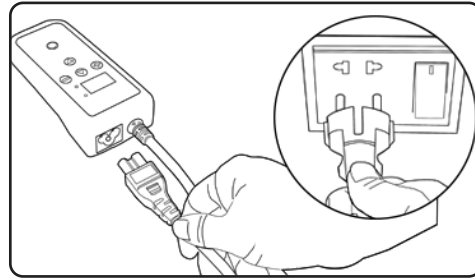

STEP 4

Press the power button on the temperature controller and set desired temperature or heating time. Wait about 10-15 minutes* and the stones will be ready to use.

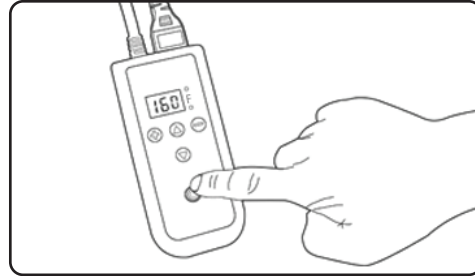

**Heating time depends on the ambient temperature and usage.*

INDICATOR STATEMENT

Status	Mode
OFF	Controller is off
ERO	NTC disconnect
Red light on	Heating
Red light twinkles	Reached the set temperature
Green light on	Countdown function is working

MAINTENANCE

1. Use a dry towel to wipe the surface of the device.
2. Place the product in the box and keep it dry if long time without using.
In order to improve user experience, there will be changes in product features without prior written notice, please refer to actual contents in the package.

WARRANTY

In order to let you use this product with confidence and better protect your rights and interests. Our company will implement after-sales service on the basis of abiding by the general administration of quality supervision, inspection and quarantine and the ministry of information industry.

The following circumstances not cover by warranty

1. The effective service or free maintenance period has been exceeded.
2. Failure to follow the correct instructions for use, maintenance, and storage leading to malfunction or damage.
3. Consumers disassemble the product themselves.
4. Product defects caused by human factors (such as operation errors, and input of inappropriate voltage and touch the water).

CAUTIONS

1. Risk of electric shock, can cause injury or death if improper operation.
2. Do not block any ventilation when using, it will cause the fire in extreme cases
3. Ensure to use in a dry environment. Do not immerse this product in water or other liquids to avoid short circuit or electric shock. Heating system is not waterproof.
4. Do not allow children or incapacitated persons to use this product.
Please keep this product out of the reach of children.
5. The Kit should be stored in a clean and dry environment.
6. Do not leave unattended.
7. This product requires professionals use. Please do not separate, maintain, or refit the product if you were not professional.
8. Heating speed may vary from the room temperature, using habits, etc.
Please set the temperature as desired avoid burns caused by excessive heat.
9. If the product has some abnormal conditions, such as bad smell, overheat, abnormal color, distortion, etc. please stop using immediately and contact the local dealer.

Features:

- Innovative Waterless Heating Technology
- Digital Thermometer Controller
- Mess-free Stone Warming, Environmentally Friendly
- With a Lightweight Energy Efficient
- Deeply Relaxing Therapeutic Experience
- Simple Electric Plug-in Design
- Quick Setup: All-in-One Portable Kit
- Portable Style Travel Case, Lightweight & Easy to Carry
- Energy Efficient Warming for Quick Heating
- Wired Temperature Controller Unit
- Configurable Time & Temp Settings
- Melt Away Stress to Soothe Muscle Tension, Soreness, Stiffness
- Heated Stone Massage Relieves Pains Caused by Sprains & Strains
- Ideal for Athletes, Home Care and Therapeutic Spa Centers

What's in the Box:

- Heating Bag
- Digital Controller
- USA Power Cable
- (8) Bamboo Massage Sticks

Technical Specs:

- Construction Material: fabric (Bag)
- Heating Time: 20-30mins
- Bamboo Maximum Heat Retention: 20-30mins
- Power Cord Length: 5.9' ft.
- Power: 110V-200W
- Rated Frequency: 60Hz
- Controller Range: 20-70°C (70-160°F)
- Total Bag Size (L x W x H): 14.6" x 10.2" x 3.5" –inches
- Bamboo Massage Sticks Dimensions: (1) 22.05" x 1.19" –inches,
 - (1) 14.96" x 1.19" –inches
 - (2) 14.96" x 0.79" –inches
 - (1) 14.96" x 0.59" -inches
 - (2) 6.3" x 1.19" -inches
 - (1) 6.3" x 1.19" -inches

