

SereneLife



VISIT US ONLINE:

Have a question?

Need service or repair?

Want to leave a comment?

PyleUSA.com/ContactUs

PYLE



SLMSGCH20

Portable Foldable Massage Chair
Programmable Heated Electric Neck
and Back Seat Massager

WARNING

To reduce the risk of burns, fire, electric shock or injury to persons:

1. The **SLMSGCH20** should never be left unattended when plugged in.
Unplug from the outlet when it is not in use and before putting on or taking off parts.
2. Close supervision is necessary when this appliance is used by or near children or disabled persons.
3. Use this appliance only for its intended use as described in this manual.
Do not use attachments not recommended by the manufacturer.
4. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
5. Keep the cord away from heated surfaces.
6. Never operate the appliance with the air openings blocked.
Keep the air openings free of lint, hair, and the like.
7. Never drop or insert any object into any opening.
8. Do not use outdoors.
9. Do not operate where aerosol spray products are being used or where oxygen is being administered.
10. To disconnect, turn the power **OFF**, then remove the plug from the outlet.
11. Do not stand on the appliance. Use only while seated.
12. The product is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.

WARNING

13. Children should be supervised to ensure that they do not play with the product.
14. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation.
The unattended use of heat by children or disabled persons may be dangerous.

FOR ADULT USE ONLY


Remove all scarves, neckties, necklaces, and jewelry before using this massager.

IMPORTANT: Any individual who may be pregnant, has a pacemaker, suffers from diabetes, phlebitis and/or thrombosis, is at an increased risk of developing blood clots, or who has pins/screws/artificial joints or other medical devices implanted in his/her body should consult with a physician before using a massaging device designed for home use.

Temperatures sufficiently high to cause burns may occur regardless of the control setting.


- Do not use on an infant or invalid or on a sleeping or unconscious person.
- Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
- Never remove cover.

LOCATION OF CONTROLS

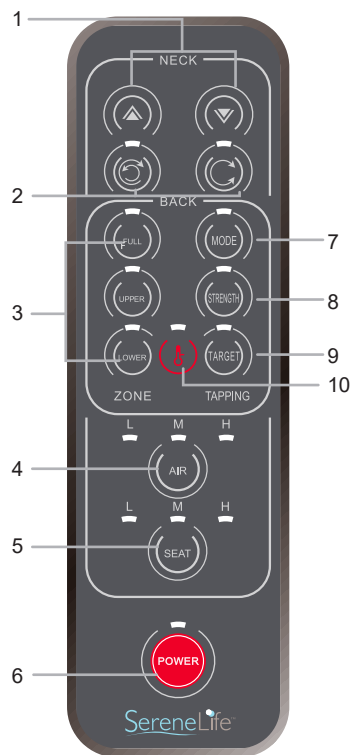

- 1. Neck Massage Adjustment**
Adjust your desired neck massage level.
- 2. Neck Massage Program**
Manual program (left button) and auto program (right button).
- 3. Back Zone Choice**
Choose your desired back massage area (FULL, UPPER, LOWER).
- 4. Airbags**
Air compression with 3 intensities.
- 5. Vibration Massage**
Seat vibration with 3 intensifies.
- 6. Power ON/OFF**
Switch ON or OFF the machine.
- 7. Modes**
Three massage modes are optional
- 8. Strength**
Three kinds of massage strength optional.
- 9. Target**
Target specific back area you like.
- 10. Heat Therapy**
Improve blood circulation.

MESSAGE UNIT

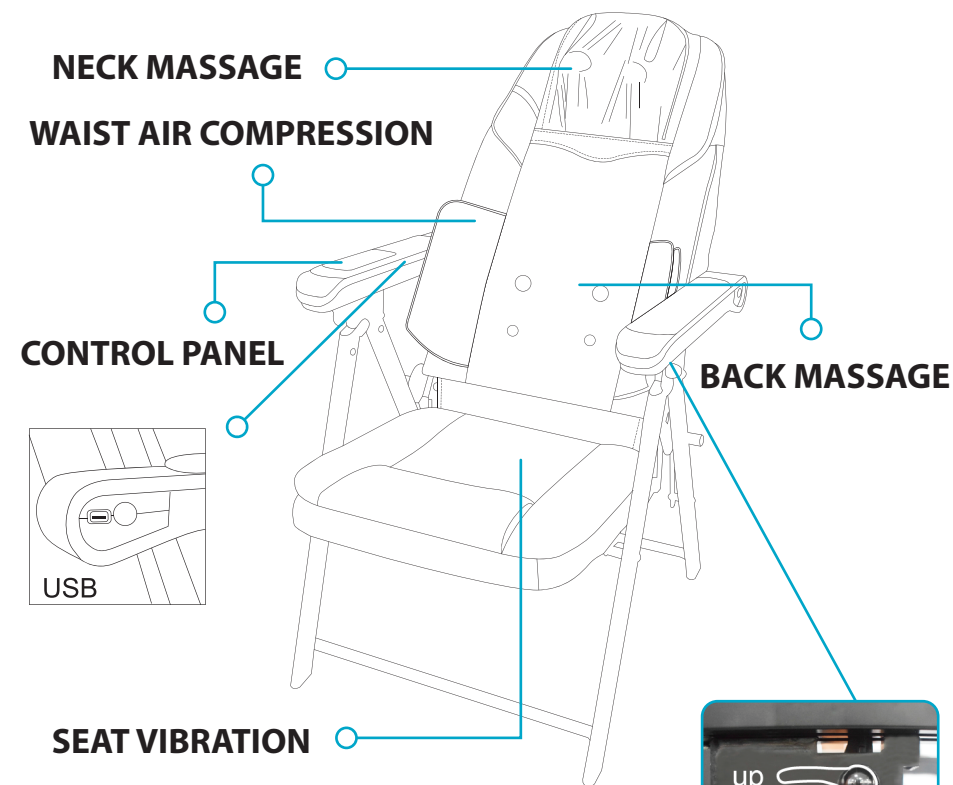

LOCK PIN INSTRUCTION

To lock pin on both sides, turn the lock pins upward. To adjust the angle or to fold the chair, turn the lock pins downward.



PREPARATION BEFORE USE POWER ON/OFF

1. Place the chair on the floor. Adjust the backrest to your desired position. Plug IN the adaptor.
2. Insert power adapter into power cord. Push until '**click**' sound is heard.
3. Insert power adapter plug into AC 120V 60HZ electric outlet.
4. Press ON the remote to turn the unit on.
5. Press again to turn the unit off.

OVERALL MASSAGE

1. Select and combine massage programs from NECK, BACK and SEAT according to your preference.
2. The selected massage programs begin.
3. During operation, you can adjust the back position, neck position, heat...etc. depending on the selected massage.
4. After the massage is completed, the rollers will return to the storage position.
5. Unplug its power plug from the electrical outlet when unit is not in use.

Note:

The unit shuts off automatically after approximately 15 minutes.

CARE AND MAINTENANCE

- Wipe dust and dirt off the Massager with a dry, soft cloth after use.
- Keep the Massager clean and away from places of high humidity.

SAFETY

1. Never leave the unit unattended, especially if children are present.
2. Never cover the unit when it is in operation.
3. Do not use this unit for more than 15 minutes at a time.
4. Extensive use at any time could lead to excessive heating that shortens the life of this unit. Should this occur, discontinue use and allow the unit to cool before operating.
5. Never use this unit directly on swollen/inflamed areas of the skin or on areas of skin eruptions. Do not use in the presence of unexplained back pain.
6. This is a non-professional unit designed to provide soothing massage to muscles. Do not use as a substitute for medical attention.
7. Do not use while in bed or before going to sleep.
The massage can create a stimulating effect and delay sleep.
8. This unit should never be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
9. Extreme heat or cold may damage the unit.
10. Place the unit in a safe, cool, dry place when not in use.
11. Do not submerge the unit in liquid.

CAUTION

- TO REDUCE THE RISK OF ELECTRIC SHOCK, DO NOT REMOVE COVER. THERE ARE NO SERVICEABLE PARTS INSIDE.
- TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THIS UNIT TO RAIN OR MOISTURE.



The lightning flash with arrow-head symbol within an equilateral triangle is intended to alert the user to the presence of uninsulated "dangerous voltage" within the unit's enclosure that may be of sufficient magnitude to constitute a risk of electric shock.



The exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the unit.

- It must only be supplied at safety extra low voltage corresponding to the marking on the appliance
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

IMPORTANT SAFETY INSTRUCTIONS

When using the massage seat topper, basic precautions should always be followed, including the following:

Read all instructions before using this massage the **SLMSGCH10**:

DANGER

To reduce the risk of electric shock:

1. Always unplug the appliance from the electrical outlet immediately after use and before cleaning.
2. Never use pins or other metallic fasteners with this appliance.
3. Keep dry — Do not operate in a water or moist condition.

TROUBLESHOOTING

Problem	Possible Cause/Action
Frictional sound of rollers against backrest cover. Motor sound during operation. Sound of air.	These sounds are the result of mechanical structure and not malfunction. This may occur during normal operation.
The Massager does not start operation.	The Massage does not start operation. Ensure that the Power is ON. Ensure that the power plug is firmly connected to the electrical outlet. Ensure power adapter is firmly connected to power cord Check that the voltage at the electrical outlet is adequate. After prolonged use, Massage motors may need to cool down to protect them.
The massage rollers stop during the massage session.	The roller mechanism may be stalled when subjected to excessive force. Release excessive force on the rollers.

SLMSGCH20

Portable Foldable Massage Chair

Programmable Heated Electric Neck and Back Seat Massager

Features:

- Relaxing Electric Massage Chair
- Space Saver Foldable Design for Easy Storage
- Selectable Massage Modes and Zones
- Vibrating Seats with Variable Intensity Level
- Head and Waist Airbags for Compression Massage
- Programmable Neck and Back Massage
- Heated Therapy Improves Blood Circulation
- Ergonomically Designed Arm Rest
- Hassle Free OSD Control Panel on the Right Arm Rest
- Enjoy a Full Back and Neck Massage!
- USB Charging Dock for Auxiliary Devices
- Power Saving: Auto Shut-Off after 15 min
- Enjoy a Full Back and Neck Massage!

What's in the Box:

- Massage Chair
- Adaptor
- User Manual

Technical Specifications:

- Construction Material: Steel, PU leather
- Maximum Support Weight: 264.5 lbs.
- Heated Seat Temperature: >86°F at 4min; 104-140°F at 15min
- Power Cord + Adaptor Line Length: 9.84" Ft.
- Power Supply: 12V 5A
- Dimension (L x W x H):
29.9" x 24.8" x 42.12" (Folded: 34.6" x 25.2" x 11.02") – inches
- Sold as: 1